



---

# Antibody Expression Vector

**pFB-CLlg-mk: Mouse Igk light chain Mammalian Expression Vector**

**Catalog#: AEV-19**

## Product Overview

This vector is a constitutive mammalian expression vector designed to deliver exceptionally high levels of antibody expression. This circular vector features an enhanced, full-length CMV promoter and other expression elements that typically enable higher expression levels. It can be used in suspension-adapted cells, such as Expi293F™ and ExpiCHO™, for transient protein expression. Additionally, it can serve as a Geneticin®-selectable expression plasmid for engineering stable cell lines. The vector carries an ampicillin resistance gene.

## Specifications

|                                  |  |
|----------------------------------|--|
| Antibiotic Resistance            | Ampicillin (AmpR)  |
| Constitutive or Inducible System | Constitutive   |
| Delivery Type                    | Transfection   |
| Promoter                         | CMV  |
| Product Type                     | Mammalian Expression Vector                                  |
| Cloning Method                   | Restriction Enzyme (5'-Agel; 3'-XhoI) or Homologous Assembly |

## Contents & Storage

- 5.0 µg of pFB-CLlg-mk in Tris-EDTA buffer
- Store at -20°C. Vectors are guaranteed stable for 12 months when properly stored.

## Materials required for antibody generation and isotype switching


- Any one of pFB-CHlg-mG1, pFB-CHlg-mG2a, pFB-CHlg-mG2b, pFB-CHlg-mG3, pFB-CHlg-mA, pFB-CHlg-mD, pFB-CHlg-mE, pFB-CHlg-mM (monomer), pFB-CHlg-mM (pentamer) expressing the mouse heavy chain.

## Vector usage



- **Entire mouse IgG, IgA, IgD, IgE, IgM monoclonal antibodies production:** cloning your heavy/light chain variable regions in the cloning sites to preserve the integrity of the heavy/light chain constant region;
- **Isotype switching:** pFB-CLlg-mk with any one of pFB-CHlg-mG1, pFB-CHlg-mG2a, pFB-CHlg-mG2b, pFB-CHlg-mG3, pFB-CHlg-mA, pFB-CHlg-mD, pFB-CHlg-mE, pFB-CHlg-mM (monomer), pFB-CHlg-mM (pentamer) expressing the mouse heavy chain could be switched a monoclonal antibody from one isotype to another, therefore, maintaining the generated new antibodies with the same antigen affinity (epitopes) but with different effector functions (enhanced or reduced or even disabled ADCC and CDC).

### Vector map

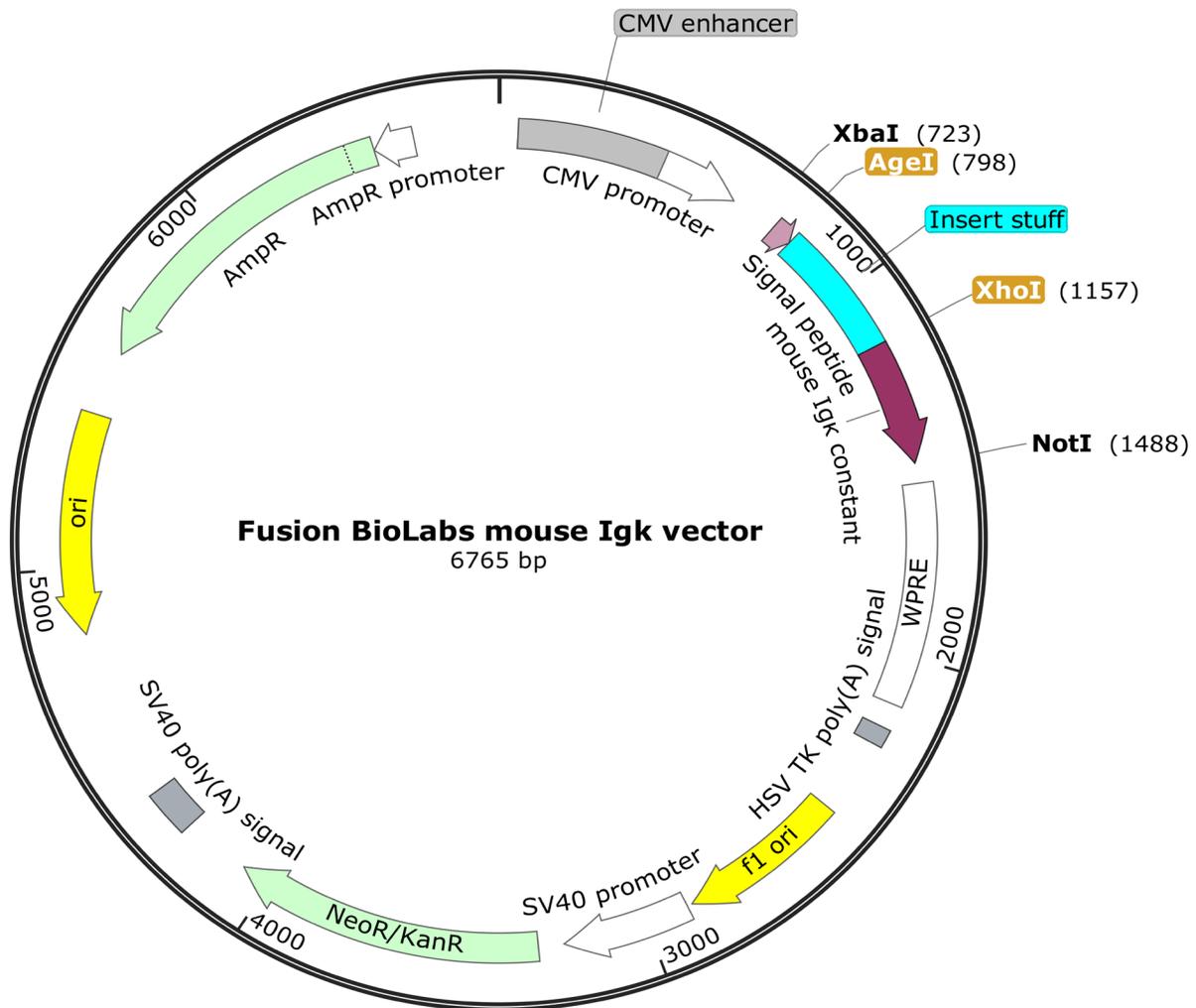



## Cloning Map



### Primer Design for Restriction Enzyme Cloning (Traditional)

Forward primer: 5'- CTAGTAGCAACTGCA**ACCGGT**GTACATTCAN<sub>VLF(15-18)</sub>

Reverse Primer: 5'- TGGTGCAGCATCAGCCCG**CTCGAG**N<sub>VLR(15-18)</sub>

**N<sub>VLF(15-18)</sub>**: append 15-18 bases encoding 5-6 amino acids starting from the beginning of the VL region to the end of the forward primer.

**N<sub>VLR(15-18)</sub>**: append the **Reverse Complement DNA Sequence** of 15-18 base encoding 5-6 amino acids of starting from the end of the VL region to the reverse primer (**LE** diamino acid will be inserted between VL and CL).

### Primer Design for Homologous Recombination Cloning (Seamless)

Forward primer: 5'- CTAGTAGCAACTGCA**ACCGGT**GTACATTCAN<sub>VLF(15-18)</sub> -3'

Reverse Primer: 5'- TGGTGCAGCATCAGCCCG**N<sub>VLR(15-18)</sub>** -3'



**N<sub>VLF</sub>(15-18):** append 15-18 bases encoding 5-6 amino acids starting from the beginning of the VL region to the end of the forward primer.

**N<sub>VLR</sub>(15-18):** append the [Reverse Complement DNA Sequence](#) of 15-18 bases encoding 5-6 amino acids starting from the end of the VL region to the reverse primer.

**Sequencing primer for confirming the correct VL insert (included in the Kit)**

pFB forward: 5'- ATGGGCGGTAGGCGTGTA-3'